



MADE IN JAPAN

Camelian *Support*

Catalogue No.	Proximal (F)	Distal (F)	Length (cm)	Inner Diameter Proximal / Distal (inch)	Tip	Hydrophilic Coating (cm)	Guidewire (inch)	Guide Catheter (inch)
CARNELIAN SUPPORT 14								
MCP090-E	2.6	1.8	90	0.021 / 0.017	Straight	60	0.014	0.035
MCP135-E	2.6	1.8	135	0.021 / 0.017	Straight	100	0.014	0.035
MCP150-E	2.6	1.8	150	0.021 / 0.017	Straight	120	0.014	0.035
CARNELIAN SUPPORT HD								
MCPH090-E	2.6	1.8	90	0.021 / 0.017	Straight	60	0.014	0.035
MCPH135-E	2.6	1.8	135	0.021 / 0.017	Straight	100	0.014	0.035
MCPH150-E	2.6	1.8	150	0.021 / 0.017	Straight	120	0.014	0.035
CARNELIAN SUPPORT DP								
MCPMS060-E	2.6	1.8	60	0.021 / 0.017	Straight	40	0.014	0.035
CARNELIAN SUPPORT BTA								
MCPB150-E	2.6	1.6	150	0.021 / 0.0155	Straight	120	0.014	0.035
CARNELIAN SUPPORT 18								
MEP070-E	2.6	1.8	70	0.021 / 0.0195	Straight	40	0.018	0.035
MEP090-E	2.6	1.8	90	0.021 / 0.0195	Straight	60	0.018	0.035
MEP135-E	2.6	1.8	135	0.021 / 0.0195	Straight	100	0.018	0.035
MEP150-E	2.6	1.8	150	0.021 / 0.0195	Straight	120	0.018	0.035

FOR ANY CASE IN PERIPHERAL INTERVENTION

Camelian *Support*



Tokai Support Catheter

Products listed may not be available in all markets. Please refer to the insert for complete instructions, contraindications, warnings and precautions.

TMP-1807284-01

Camelion Support

Multistep Shaft with Optimal Balance

provides additional strength, pushability and crossability

Tungsten Braid and Platinum Marker Band

provide excellent visibility enabling precise positioning

Hydrophilic Coating


facilitates easy and smooth tracking

Extensive Selection

available for any case in peripheral intervention

The Proven Quality and Performance


the most accepted support catheter in Japan

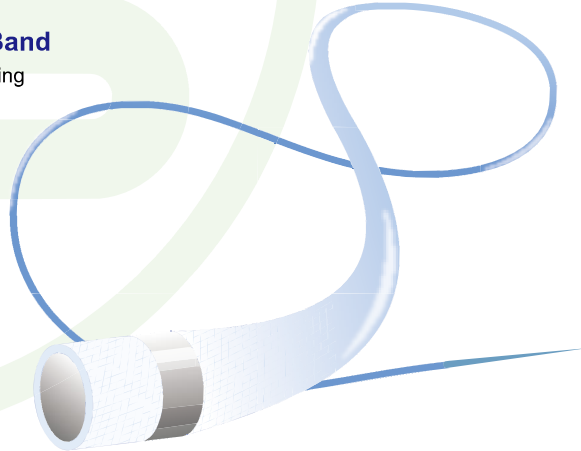

Recommended Indications

	14	HD	DP	BTA	18
Antegrade approach	✓	✓			✓
Retrograde approach	✓	✓			✓
Trans-collateral approach					
SFA-CTO		✓			✓
Below the Knee				✓	
Distal puncture			✓		

Stiffness

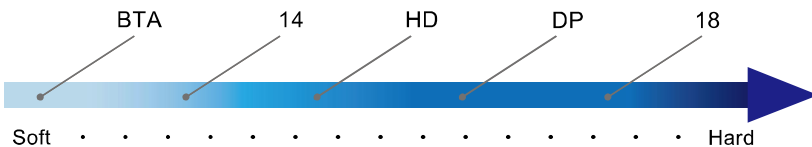


14 Support 14

Workhorse

1.8Fr/2.6Fr 90, 135, 150cm



- Well-balanced and Multi-purposed support catheter

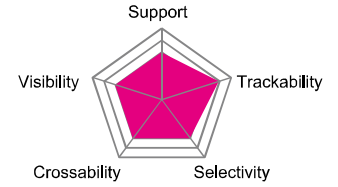


HD Support HD

Hard shaft

1.8Fr/2.6Fr 90, 135, 150cm



- Hard shaft provides additional back-up force

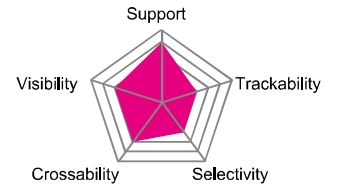


DP Support DP

Distal Puncture

1.8Fr/2.6Fr 60cm



- Specially designed for Distal Puncture
- Short tapered tip provides easy insertion and prevents blood from leaking at puncture site

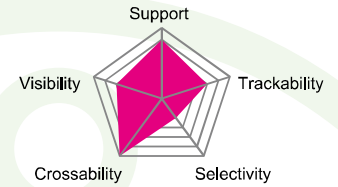


BTA Support BTA

Below The Ankle

1.6Fr/2.6Fr 150cm



- 300mm of 1.6Fr. flat shaft allows tracking through long and tortuous anatomy such as trans-collateral approach

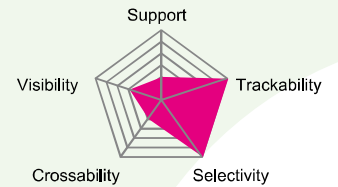

18 Support 18

0.018" compatible

1.8Fr/2.6Fr 70, 90, 135, 150cm



- Low profile tapered tip with additional strength of shaft facilitates easy crossing of challenging lesion.

